

**FORTRESS**

Nothing is Built like a Fortress

# LUXURY IN MOTION

**S700**



# S700

## PERFORMANCE LUXURY SCOOTER

---

The Fortress S700, with superior performance, features the perfect companion for all your adventures.

There's no need to worry about distance either as the large battery capacity will power an extended range and a lively top speed of **15 km/h**.

The high performance design means varied terrains are tackled with ease, and the large tires and advanced 4-wheel suspension provides a smooth, comfortable ride each and every time.







Max Speed: **15 km/h**



Colour: **Carbon Metallic Grey**



Max. Battery Capacity: **75 Amps (GP24)**



Wheel Size: **33 cm (13 inches)**



Max. User Weight: **159 kg. (350 lbs.)**

# DELUXE LED LIGHT PACKAGE

---

**Powerful all-around LED lighting. With lifetime performance, changing bulbs is a thing of the past.**

# INNOVATION IN DESIGN

---

With advanced LED lighting technology, the S700 scooter is 400 times more energy efficient than standard bulbs, meaning more battery power and range every trip.

**Safety is standard! S700 features high visibility LED turn signals, 4 way flashers and automatic brake lights.**






**Innovative active anti-tip wheels**  
move up as you ride over obstacles, and  
keep you safe on slopes.

# LUXURY SEATING

---

**The multi-adjustable luxury seating  
is designed to make every journey a  
comfortable one.**

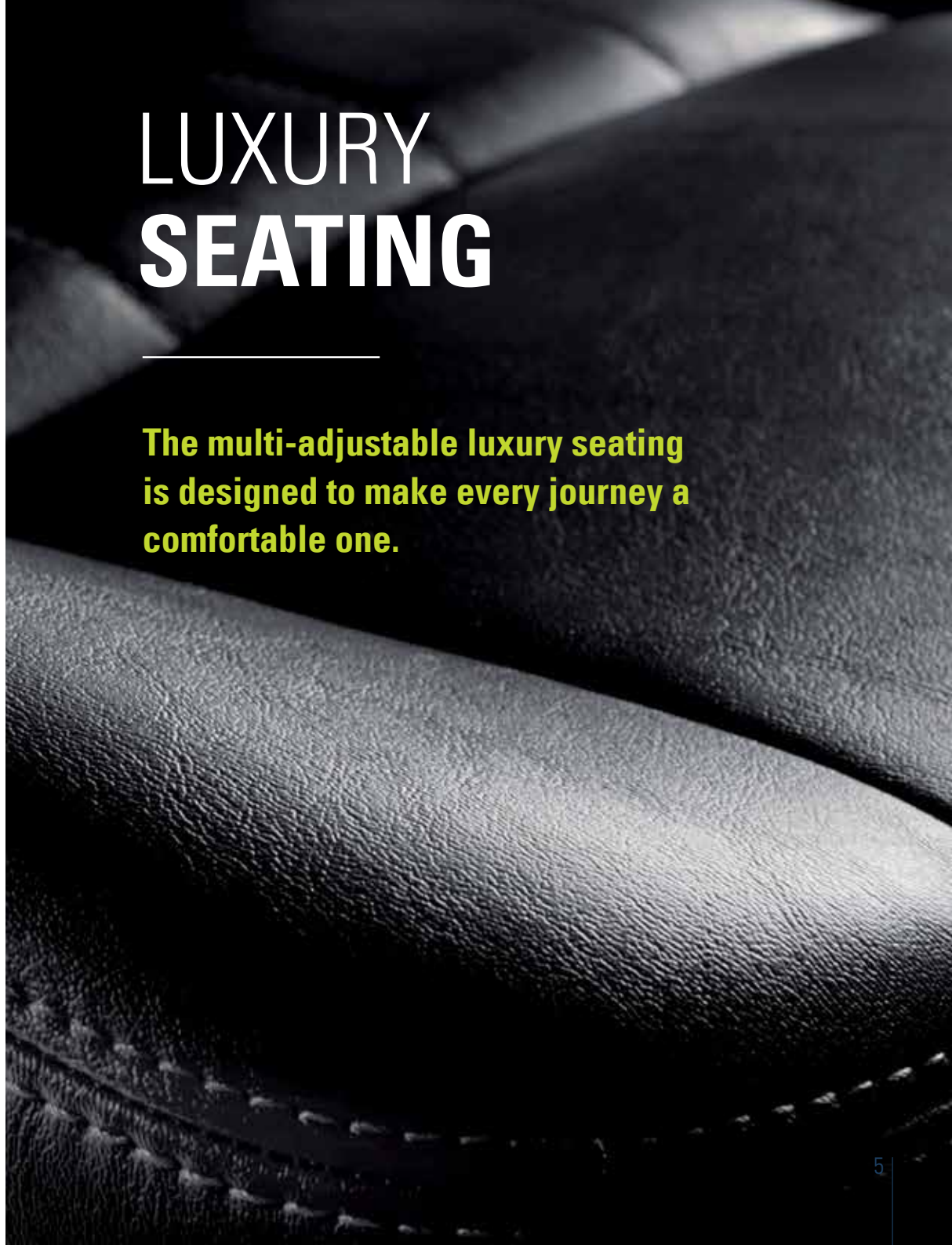

# COMFORT ON EVERY TRIP

---

The S700 advanced seating solution takes comfort to the highest level. From the seat height, seat depth and back recline adjustments, to the flip-up armrests including width, angle and depth adjustments – every feature has been designed for comfort – for you.



→ Infinite seat depth adjustment



→ Quick backrest angle adjustment



→ Seat rotation enables comfortable transfers



# CONTROLS DESIGNED FOR YOU

Finding your perfect driving position is a breeze with the infinitely adjustable tiller. Clever design means the headlight automatically tilts as the tiller is adjusted.

Soft touch controls featuring easy to use speed adjustment, LED battery gauge and a powerful horn.



Unique forward/reverse lever design means natural driving using finger or thumb control. The large ergonomically designed delta handlebar makes every journey a comfortable one.

Safety and style are standard. Every S700 scooter comes with two carbon effect mirrors.



Deluxe extra large basket comes standard to ensure plenty of room for your belongings.

S700



# TECHNICAL DATA

	<i>S700</i>
<b>Max Speed</b>	15 km/h (9.3 mph)
<b>Max. User Weight</b>	159 kg. (350 lbs.)
<b>Max Battery Capacity</b>	75 Amps (GP24)
<b>Charger</b>	8 A off-board
<b>Max. Potential Range*</b>	55 km (34 miles)
<b>Max. Power (peak)</b>	2700 W (4-pole motor)
<b>Deluxe LED Light Package</b>	Yes
<b>Easy Adjust Tiller &amp; Throttle Control</b>	Yes
<b>Suspension</b>	Deluxe Front & Rear Suspension (Front Independent suspension)
<b>Tires</b>	13" pneumatic
<b>Overall Width / Length</b>	66 cm / 149 cm (26" / 59")
<b>Seat Width x Depth</b>	48 cm x 48 cm (19" x 19")
<b>Seat Height Range</b>	45 cm - 50 cm (18" - 20")
<b>Minimum Turning Distance</b>	2.2 m (7 feet)
<b>Colour</b>	Carbon Metallic Grey
<b>Scooter Weight w/ Batteries 75 Amp</b>	320 lbs. inc. GR24

\* Max potential range based on industry benchmark range figures for similar products and battery sizes



All over Canada, thousands of people have made their lives easier with a FORTRESS mobility scooter by Sunrise Medical. For over 30 years, FORTRESS has designed and built mobility scooters to the highest standards, giving the comfort and safety every customer needs. For further information on product options please check the order form on our website. You can view this document also in a larger font on a PDF file at [www.SunriseMedical.ca](http://www.SunriseMedical.ca)



## Sunrise Medical Canada Inc.

237 Romina Drive, Unit 3  
 Concord, ON  
 Tel.: 1.800.263.3390  
 Fax: 1.800.561.5834  
 E-mail: [mkt.canada@sunmed.com](mailto:mkt.canada@sunmed.com)  
[www.SunriseMedical.ca](http://www.SunriseMedical.ca)